

Is Achieving
Primary Stability
A Challenge?
Try The OSSEOTITE^{®2}
Parallel Walled
Implant!



Comprehensive
Treatment

Enhanced
Knowledge

Business
Excellence



Introducing **NEW** OSSEOTITE^{®2} Implants


The Solution Designed To Help Achieve Primary Stability With Parallel Walled Implants:

- New Design Features For A Tighter Osteotomy Fit
- Increased Lateral Threads For "Bite-In-Bone" Engagement
- Designed To Provide More Surface Area In Direct Contact With Bone



Providing Clinicians One Solution At A Time With OSSEOTITE^{®2} Implants

- More Surface Area Than The Previous Design For Greater Immediate Bone-To-Implant Contact
- Implant Design Allows For Steady Increase Of Insertional Torque Throughout Placement Indicating Mechanical Engagement Into Bone

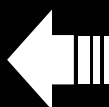

**Scan With Your
Smartphone!**

In order to scan QR codes,
your mobile device
must have a QR code
reader installed.

BIOMET 3i[™]
PROVIDING SOLUTIONS - ONE PATIENT AT A TIME[™]


Subcrestal Placement For OSSEOTITE² Certain[®] PREVAIL[®], OSSEOTITE² Certain[®] And OSSEOTITE² External Hex Parallel Walled Implants

NOTE:


- The recommended drill speed for drills 3.85mm diameter or less is 1200–1500rpm.
- The recommended drill speed for drills 4.25mm diameter or greater is 900rpm.
- The implant placement torque may exceed 50Ncm.
- The recommended implant placement speed is 15–20rpm.
- Final Twist Drill selection is based on clinician evaluation of bone quality.
- Hand ratcheting may be necessary to fully seat the implant into the osteotomy.
- **Tapping is required in dense (Type I) bone for 5.0mm, 6.0mm, 5/4mm and 6/5mm diameter implants.**

D = Diameter
C = Collar
P = Platform
L = Length

ACT[®] Twist Drill Depth Marks

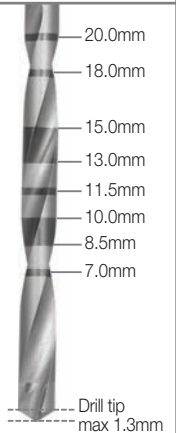


OSSEOTITE² Certain[®] 3.25mm And OSSEOTITE² External Hex 3.25mm Implants

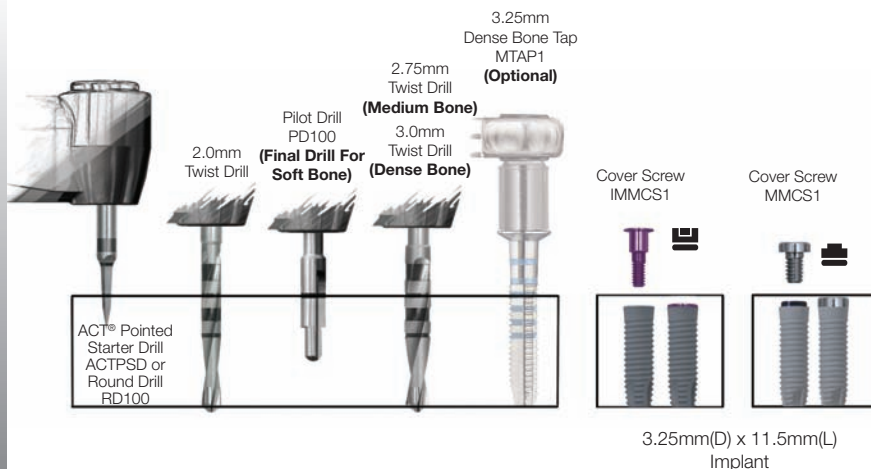


IMPORTANT NOTE:

Exceeding an insertion torque of more than 90Ncm may deform or strip the driver tip or the implant's internal hex and may possibly delay the surgical procedure.


OSSEOTITE² External Hex 3.75mm Implant

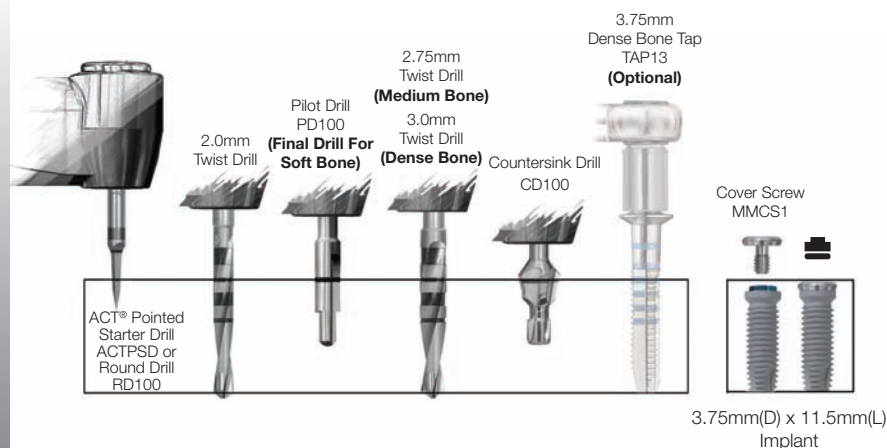


OSSEOTITE² Certain[®] PREVAIL[®] 4/3mm, OSSEOTITE² Certain[®] 4.0mm And OSSEOTITE² External Hex 4.0mm Implants

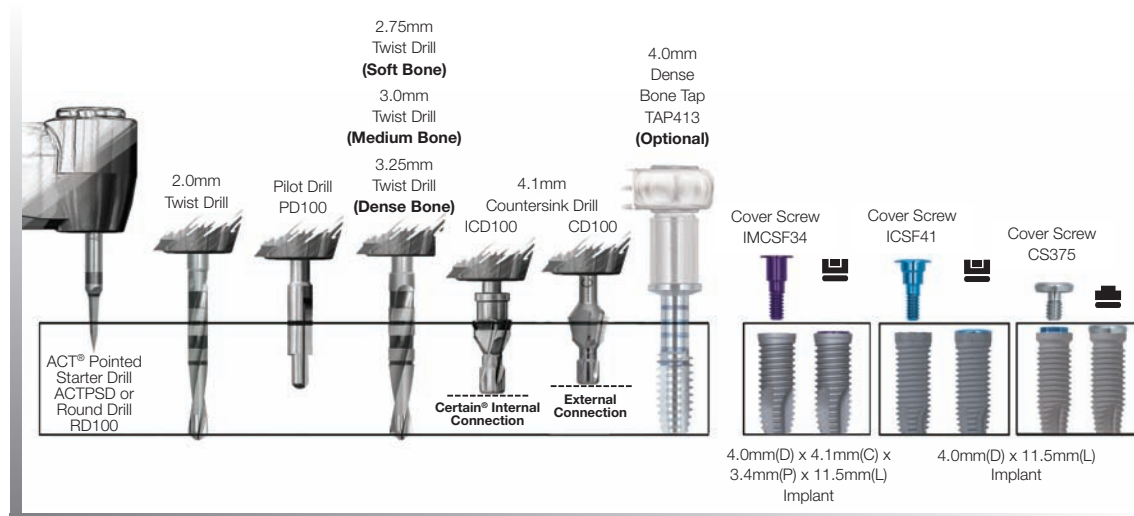

OSSEOTITE² Certain[®] PREVAIL[®] 5/4mm, OSSEOTITE² Certain[®] 5.0mm And OSSEOTITE² External Hex 5.0mm Implants

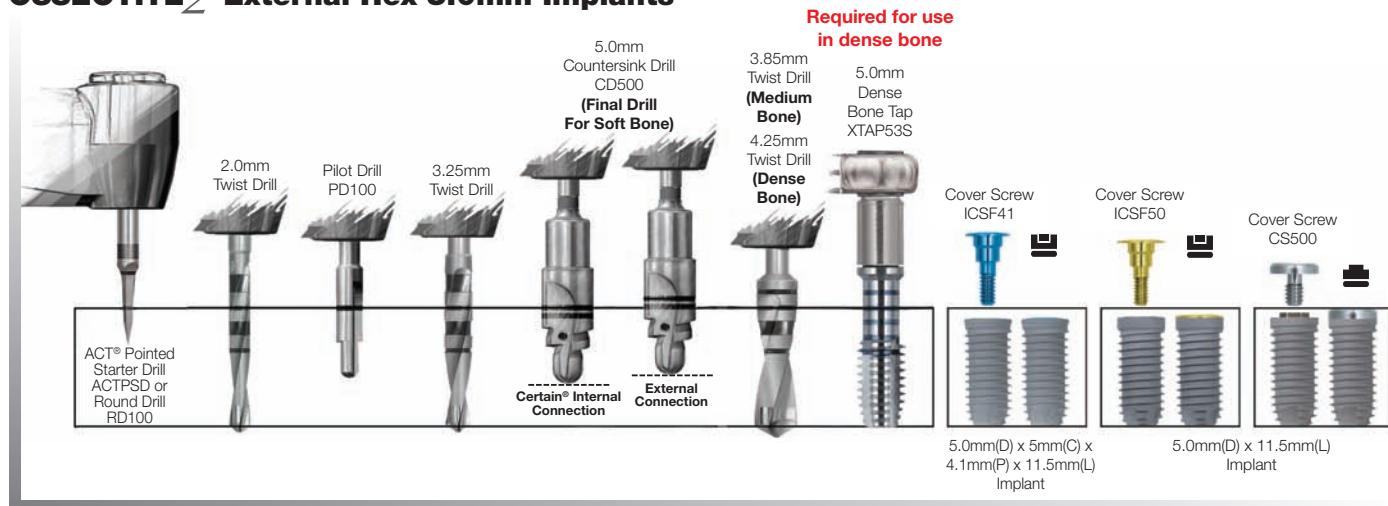

OSSEOTITE² Certain[®] PREVAIL[®] 6/4mm, OSSEOTITE² Certain[®] 6.0mm And OSSEOTITE² External Hex 6.0mm Implants

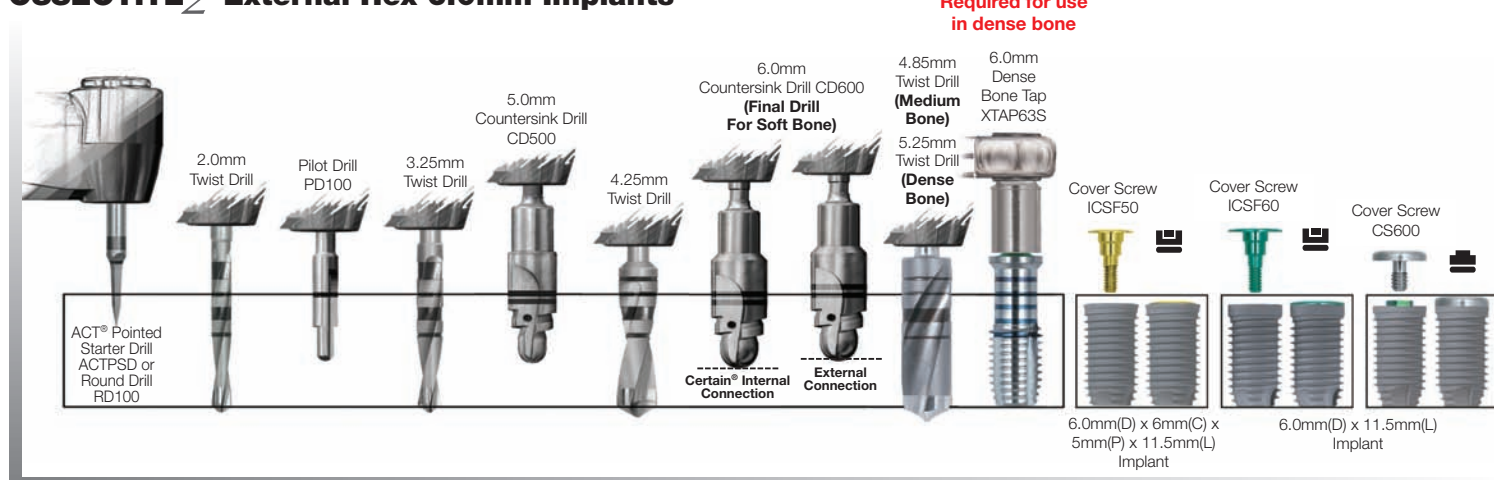

To Find Out How You Can Achieve Primary Stability With Parallel Walled Implants, Please Contact Your Local **BIOMET 3i** Sales Representative Today

SUBSIDIARIES

AUSTRALIA

Phone: +61-2-9855-4444
Fax: +61-2-9888-9900

AUSTRIA

Phone: 0800070017
Fax: 0800070018

BELGIUM

Phone: 0800 50311
Fax: 0800 50312

BRAZIL

Phone: +55-11-3568-1300
Fax: +55-11-5081-7484

CANADA

Phone: +514-956-9843
Fax: +514-956-9844

DENMARK

Phone: 0045 33 12 70 08
Fax: 0045 33 12 70 03

FRANCE

Phone: 01 41 05 43 43
Fax: 01 41 05 43 40

GERMANY

Phone: 0800 1016420
Fax: 0800 3131111

IRELAND

Phone: 1800 552 752
Fax: 1800 656 608

JAPAN

Phone: +81-66868-3012
Fax: +81-66868-2444

KOREA

Phone: +82-2-5678-550
Fax: +82-2-5678-577

MEXICO

Phone: +52-55-2282-0120
Fax: +52-55-2282-0120 ext. 20

THE NETHERLANDS

Phone: 078 62 92 800
Fax: 078 62 92 801

NEW ZEALAND

Phone: +64-508-122-221
Fax: +64-508-133-331

NORDIC REGION HEADQUARTERS

Phone: +46-40-17-6090
Fax: +46-40-17-6099

NORWAY

Phone: 0047 22 20 12 70
Fax: 0047 22 83 97 40

PORTUGAL

Phone: +34-902-34-34-31
Fax: +34-93-445-81-36
+351-21-000-16-75

SPAIN

Phone: +34-902-34-34-31
Fax: +34-93-445-81-36

SWEDEN

Phone: 0046 (0)20 0125509
Fax: 0046 (0)20 0125506

SWITZERLAND

Phone: 0800 246638
Fax: 0800 246639

U.K.

Phone: 0800 652 1233
Fax: 0800 840 6814

DISTRIBUTORS

ARGENTINA

Dentalmax SA
Phone: +541-1482-71001
Fax: +541-1482-67373

CHILE

Cybel SA
Phone: +56-2-2321883
Fax: +56-2-2330176

CHINA

Atek Inc.
Phone: +86-21-6329-1265
Fax: +86-21-6329-1620

COLOMBIA

3i Colombia
Phone: +571-612-9362
Fax: +571-620-6412

COSTA RICA

Implantec S.A.
Phone: +506-234-9043
Fax: +506-224-7620

EL SALVADOR

Dentimer SA de CV
Phone: +503-263-6350
Fax: +503-263-6676

GREECE

Impladent Dental Implants, LLC
Phone: +30-2310-501-651
Fax: +30-2310-862-090

ISRAEL

H.A. Systems
Phone: +972-3-6138777
Fax: +972-3-6138778

HONG KONG

Ositek Inc., Ltd.
Phone: +852-8121-6601
Fax: +852-3747-3754

ITALY

Biomax spa
Phone: +39-0444-913410
Fax: +39-0444-913695

MIDDLE EAST

3i MENA s.a.l.
Middle East And North Africa
Phone: +961-1-694000
Fax: +961-1-694222

PARAGUAY

Andres H. Arce y Cia SRL
Phone: +595-21-208185
Fax: +595-21-496291

POLAND

Dental Depot Wasio
Phone: +48 71 335 70 71
Fax: +48 71 335 70 90

RUSSIA

Com-Dental
Phone: +7-495-580-3080
Fax: +7-495-580-3081

SINGAPORE

Asia Implant Support & Services
Phone: +65-6223-2229
Fax: +65-6220-3538

SOUTH AFRICA

Selective Aesthetics
Phone: +27-11-991-7068
Fax: +27-11-672-1391

TAIWAN

Kuo Hwa Dental Suppliers Co., Ltd.
Phone: +886-2-2226-1770
Fax: +886-2-2226-8747

THAILAND

3i (Thailand) Co., LTD.
Phone: +662-252-6685
Fax: +662-252-6686

TURKEY

Tamer Med. A.Ş.-Turkey
Phone: +90-212-465-3352
Fax: +90-212-465-3502

UKRAINE

Victoria Premium
Phone: +38-044-331-10-09
Fax: +38-044-331-10-70

URUGUAY

Pro3implant S.R.L.
Phone: +598-2-4034163
Fax: +598-2-4034163

VIETNAM

D.O.E.
(Dentistry Of Excellence Vietnam)
Phone: +84 89 25 37 03
Fax: +84 89 25 37 00

ACT, Certain, OSSEOTITE and PREVEAL are registered trademarks and Providing Solutions – One Patient At A Time and design are trademarks of BIOMET 3i LLC. BIOMET 3i and design are trademarks and BIOMET is a registered trademark of BIOMET, Inc. ©2011 BIOMET 3i LLC. All rights reserved.



REV A 05/11

To Receive Information About BIOMET 3i's Products, Services And Events By Email, Please Visit www.biomet3i.com

Not Available In All Markets, Please Consult Your Local BIOMET 3i Sales Representative For Availability

BIOMET 3i

4555 Riverside Drive
Palm Beach Gardens, FL 33410
1-800-342-5454
Outside The U.S.: +1-561-776-6700
Fax: +1-561-776-1272
www.biomet3i.com



BIOMET 3i

Europe, Middle East & Africa
WTC Almeda Park, Ed. 1, Planta 1ª
Pl. de la Pau, s/n
08940, Cornellà de Llobregat
(Barcelona) Spain
Phone: +34-93-470-55-00
Fax: +34-93-371-78-49

BIOMET 3i
PROVIDING SOLUTIONS – ONE PATIENT AT A TIME™

Follow us on:

